

University of Dundee

**Ketamine inhibits synaptic transmission and nicotinic acetylcholine receptor-mediated responses in rat intracardiac ganglia in situ**

Harper, Alexander A.; Rimmer, Katrina; Dyavanapalli, Jhansi; McArthur, Jeffrey R.; Adams, David J.

*Published in:*  
Neuropharmacology

*DOI:*  
[10.1016/j.neuropharm.2019.107932](https://doi.org/10.1016/j.neuropharm.2019.107932)

*Publication date:*  
2020

*Licence:*  
CC BY-NC-ND

*Document Version*  
Peer reviewed version

[Link to publication in Discovery Research Portal](#)

*Citation for published version (APA):*

Harper, A. A., Rimmer, K., Dyavanapalli, J., McArthur, J. R., & Adams, D. J. (2020). Ketamine inhibits synaptic transmission and nicotinic acetylcholine receptor-mediated responses in rat intracardiac ganglia *in situ*. *Neuropharmacology*, 165, [107932]. <https://doi.org/10.1016/j.neuropharm.2019.107932>

**General rights**

Copyright and moral rights for the publications made accessible in Discovery Research Portal are retained by the authors and/or other copyright owners and it is a condition of accessing publications that users recognise and abide by the legal requirements associated with these rights.

- Users may download and print one copy of any publication from Discovery Research Portal for the purpose of private study or research.
- You may not further distribute the material or use it for any profit-making activity or commercial gain.
- You may freely distribute the URL identifying the publication in the public portal.

**Take down policy**

If you believe that this document breaches copyright please contact us providing details, and we will remove access to the work immediately and investigate your claim.

## Supplementary Material

### **Ketamine inhibits synaptic transmission and nicotinic acetylcholine receptor-evoked responses in rat intracardiac ganglia *in situ***


Alexander A. Harper, Katrina Rimmer, Jhansi Dyavanapalli, Jeffrey R. McArthur and David J. Adams\*

#### **Supplementary Methods**

##### ***Electrophysiology***

HEK293T cells (ATCC, USA) were cultured in DMEM (Invitrogen Life Technologies, Australia) supplemented with 10% FBS (Bovigen, Australia), 1% penicillin and streptomycin and 1x GlutaMAX (Invitrogen Life Technologies) at 37°C in 5% CO<sub>2</sub>. HEK293T cells were plated on 12 mm glass coverslips and transiently transfected using the calcium phosphate method. Plasmid cDNAs encoding human Cav2.1 (provided by Dr J. Striessnig, University of Innsbruck, Innsbruck, Austria), Cav2.2 (provided by Dr D.T. Yue, The John Hopkins University School of Medicine, Baltimore, MD, USA) or Cav2.3 (provided by Dr T. Schneider, University of Cologne, Cologne, Germany) were co-transfected with human  $\alpha 2\delta 1$  and  $\beta 3$  together with a green fluorescent protein (GFP) expressing plasmid (for identification of transfected cells).

Whole-cell patch clamp recordings were performed within 24-48 hours post-transfection at room temperature (21-23°C). Cells were constantly superfused with extracellular solution containing (in mM): 100 NaCl, 10 CaCl<sub>2</sub>, 1 MgCl<sub>2</sub>, 5 CsCl, 30 TEA-Cl, 10 Glucose, 10 HEPES (pH 7.35 with TEA-OH; ~310 mOsmol·kg<sup>-1</sup>). Fire-polished borosilicate pipettes (1-3 M $\Omega$ ) were filled with intracellular solution containing (in mM): 140 K-Gluconate, 5 NaCl, 2 MgCl<sub>2</sub>, 5 EGTA, 10 HEPES (pH 7.2 with KOH; ~290 mOsmol·kg<sup>-1</sup>). Calcium currents were elicited by a test depolarization to +30 mV (50 ms duration) from a holding potential of -90 mV at 0.1 Hz. Use of ketamine HCl (Ketamil) was approved by NSW Ministry of Health (Approval #A-201911-74).

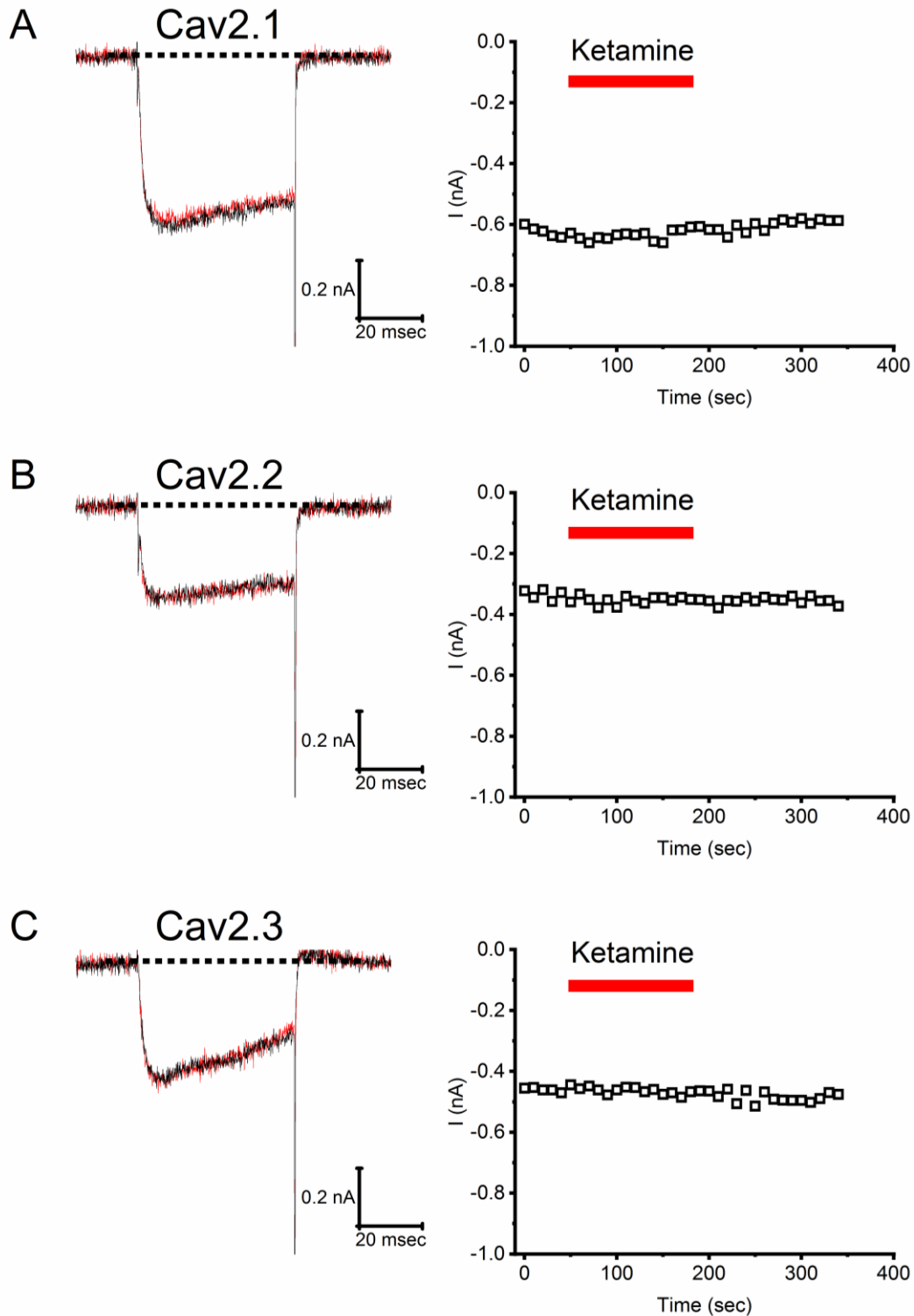

**Figure S1.** Ketamine (100  $\mu$ M) does not inhibit human Cav2.1, Cav2.2 or Cav2.3 channels heterologously expressed in HEK293T cells. Representative depolarization-activated  $\text{Ca}^{2+}$  currents (right) elicited from a holding potential of  $-90$  mV to a test potential of  $+30$  mV (50 msec duration; 0.1 Hz) from Cav2.1 (A), Cav2.2 (B) and Cav2.3 (C) in the absence (control, black) and presence of 100  $\mu$ M ketamine (red). Corresponding timeplots of  $\text{Ca}^{2+}$  current amplitude before, during and upon washout of 100  $\mu$ M ketamine application for each Cav are shown in the LHS panels.