



University of Dundee

Sperm Toolbox-A selection of small molecules to study human spermatozoa

Gruber, Franz S.; Richardson, Anthony; Johnston, Zoe C.; Myles, Rachel; Norcross, Neil R.; Day, David P.

Published in:
PLoS ONE

DOI:
[10.1371/journal.pone.0297666](https://doi.org/10.1371/journal.pone.0297666)

Publication date:
2024

Licence:
CC BY

Document Version
Publisher's PDF, also known as Version of record

[Link to publication in Discovery Research Portal](#)

Citation for published version (APA):

Gruber, F. S., Richardson, A., Johnston, Z. C., Myles, R., Norcross, N. R., Day, D. P., Georgiou, I., Sesma-Sanz, L., Wilson, C., Read, K. D., Martins da Silva, S., Barratt, C. L. R., Gilbert, I. H., & Swedlow, J. R. (2024). Sperm Toolbox-A selection of small molecules to study human spermatozoa. *PLoS ONE*, *19*(2), Article e0297666. <https://doi.org/10.1371/journal.pone.0297666>

General rights

Copyright and moral rights for the publications made accessible in Discovery Research Portal are retained by the authors and/or other copyright owners and it is a condition of accessing publications that users recognise and abide by the legal requirements associated with these rights.

Take down policy

If you believe that this document breaches copyright please contact us providing details, and we will remove access to the work immediately and investigate your claim.

RESEARCH ARTICLE

Sperm Toolbox—A selection of small molecules to study human spermatozoa

Franz S. Gruber¹ , Anthony Richardson² , Zoe C. Johnston³ , Rachel Myles³, Neil R. Norcross², David P. Day², Irene Georgiou², Laura Sesma-Sanz¹, Caroline Wilson², Kevin D. Read², Sarah Martins da Silva², Christopher L. R. Barratt³, Ian H. Gilbert², Jason R. Swedlow¹ *

1 Divisions of Computational Biology and Molecular, Cell and Developmental Biology, and National Phenotypic Screening Centre, School of Life Sciences, University of Dundee, Dundee, United Kingdom, **2** Drug Discovery Unit, Division of Biological Chemistry and Drug Discovery, School of Life Sciences, University of Dundee, Dundee, United Kingdom, **3** Division of Systems Medicine, School of Medicine, Ninewells Hospital and Medical School, University of Dundee, Dundee, United Kingdom

 These authors contributed equally to this work.

* jrswedlow@dundee.ac.uk


 OPEN ACCESS

Citation: Gruber FS, Richardson A, Johnston ZC, Myles R, Norcross NR, Day DP, et al. (2024) Sperm Toolbox—A selection of small molecules to study human spermatozoa. PLoS ONE 19(2): e0297666. <https://doi.org/10.1371/journal.pone.0297666>

Editor: Wagdy Mohamed Eldehna, Kafrelsheikh University Faculty of Pharmacy, EGYPT

Received: September 7, 2023

Accepted: January 11, 2024

Published: February 20, 2024

Peer Review History: PLOS recognizes the benefits of transparency in the peer review process; therefore, we enable the publication of all of the content of peer review and author responses alongside final, published articles. The editorial history of this article is available here: <https://doi.org/10.1371/journal.pone.0297666>

Copyright: © 2024 Gruber et al. This is an open access article distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited.

Data Availability Statement: All process data are included in the Supplemental Info. Data have also been deposited at zenodo, and are available at

Abstract

Male contraceptive options and infertility treatments are limited, and almost all innovation has been limited to updates to medically assisted reproduction protocols and methods. To accelerate the development of drugs that can either improve or inhibit fertility, we established a small molecule library as a toolbox for assay development and screening campaigns using human spermatozoa. We have profiled all compounds in the Sperm Toolbox in several automated high-throughput assays that measure stimulation or inhibition of sperm motility or the acrosome reaction. We have assayed motility under non-capacitating and capacitating conditions to distinguish between pathways operating under these different physiological states. We also assayed cell viability to ensure any effects on sperm function are specific. A key advantage of our studies is that all compounds are assayed together in the same experimental conditions, which allows quantitative comparisons of their effects in complementary functional assays. We have combined the resulting datasets to generate fingerprints of the Sperm Toolbox compounds on sperm function. The data are included in an on-line R-based app for convenient querying.

Introduction

One of the grand challenges for science and society is the development of novel contraceptives that are safe and deliver control of reproduction to men and women all over the world. In parallel, a deeper understanding of sperm maturation, motility, fertilization, and possible interventions for male infertility is required to counter declining birth rates. Recently, the application of advanced technologies has helped move the human fertility field forward, revealing new potential targets for both contraception and fertility, and new compounds as contraceptive candidates (reviewed by [1–3]). A key part of this effort is the development of mechanistic insight and drug candidates that affect the behaviour and properties of human

<https://doi.org/10.5281/zenodo.10156888> (DOI: 10.5281/zenodo.10156888).

Funding: The authors acknowledge support by the Bill & Melinda Gates Foundation (INV-007117). The funders had no role in study design, data collection and analysis, decision to publish, or preparation of the manuscript.

Competing interests: The authors have declared that no competing interests exist.

spermatozoa [4–7]. Towards this end, we have recently developed an automated live human sperm, target-agnostic, phenotypic screening platform that can profile the effects of 1000s of compounds on spermatozoa. These parameters are then used to identify modifiers of sperm function [8, 9]. These modifiers can cause either positive or negative effects and therefore are starting points for either fertility treatments (including medical assisted reproductive technologies) or contraceptive development.

A limitation of the phenotypic assays we have employed is the lack of mechanistic information on compounds that score in the assay. Because there is so little known about the molecular mechanisms that control sperm motility, we are unable to develop hypotheses of potential targets when new hits are identified in the sperm motility assay. However, it is now well-established that features derived from images in phenotypic assays, when compared to phenotypes measured on reference compounds, can accurately stratify genetic or chemical perturbations in phenotypic screens [10, 11]. We have therefore collected 83 reference compounds that have been reported in the literature to affect sperm function in various assays as well as compounds we have observed to modify sperm function in our own assays, to form the Sperm Toolbox (STB) (Table 1). In constructing the STB, we have focussed on compounds reported in the literature that affect human systems, but compounds tested in relevant animal models, cell-based, or biochemical systems have also been included (Fig 1A). Published analyses of the compounds in the STB have indicated effects on sperm function e.g. motility, capacitation, acrosome reaction, Ca²⁺ signalling, viability, or related to other processes (e.g. spermatogenesis, liquefaction, or ejaculation; Fig 1C). Compound annotations are included in the dataset derived from the STB and indicate that Enzymes/Kinases, Receptors, Channels and Transporters are the most represented target of the STB compound classes (Fig 1D). Similar toolboxes have been established for other diseases or targets, e.g. the MMV Pandemic Response Box [12], the SGC Epigenetic Chemical Probe Collection [13], or the Drug Repurposing Hub [14].

List of compounds included in Sperm Toolbox. Published Phenotypes are summaries of reported effects on sperm function related to other processes (e.g. spermatogenesis, liquefaction, or ejaculation). Compounds from our own screening campaigns, which have not been published and are included into the Sperm Toolbox have been labelled with ‘This study’ (Fig 1B).

To demonstrate the value of the STB, we have assayed STB compounds in a series of high throughput phenotypic assays of sperm function (Fig 1A). The resulting quantitative fingerprint of sperm function modifiers provides a useful reference for drug discovery campaigns as well as experiments that explore relevant molecular mechanisms. An advantage of using a high-throughput platform is the ability to simultaneously compare many compounds under the same assay conditions (same buffer system, same sample), within a short amount of time. This overcomes limitations of screening compounds at low throughput, where only few compounds can be tested with the same donor pool, one compound at a time. Invariably some of the compounds we have included in the STB have proven to have little effect in human sperm, at least under the conditions and at the concentrations we have used. We have maintained these compounds in the STB and reported the results to ensure the published literature is correct and authoritative.

Results and discussion

We tested the STB compounds in multiple functional sperm assays, as well as more general cytotoxicity and aqueous solubility assays (Fig 1A). Using information from the Drug Repurposing Hub [14], we generated a protein-protein interaction network of targets related to STB compounds (Fig 1E). This network consists of multiple subclusters related to functions in spermatozoa (Fig 1F) and other biological processes, functions, and pathways (Fig 1G).

Table 1. Sperm Toolbox compounds.

STB Number	Name	Published Phenotype	References
1	ASI-8	Blocks soluble adenylyl cyclase (sAC)	[15–17]
2	bpV(phen)	Blocks progesterone induced acrosome reaction; blocks progesterone induced Ca ²⁺ response	[18]
3	Tideglusib	Spermicidal; vaginal treatment	[19]
4	Genistein	Blocks progesterone induced acrosome reaction; blocks progesterone induced Ca ²⁺ response	[18]
5	T16Ainh-A01	Blocks acrosome reaction induced by recombinant ZP3	[20]
6	Niflumic acid	Blocks acrosome reaction induced by recombinant rhZP3	[20]
7	Hexachlorophene	Reduces sperm motility; blocks sAC	[15]
8	Fenvalerate	Blocks progesterone induced acrosome reaction	[18]
9	Niclosamide	Shape abnormality	[21]
10	BMS-189453	Interferes with spermatogenesis	[22]
11	Dipyridamole	Enhances sperm motility	[23]
12	Cloperastine hydrochloride	Reduces sperm motility	This study
13	Resiquimod	Reduction of sperm motility; only sperm with X chromosome	[8, 24]
14	Compound B4	Reduces sperm motility	[25]
15	2-Guanidinobenzimidazole	Reduces sperm motility	[26, 27]
16	Alexidine dihydrochloride	Reduces sperm motility	[8]
17	TP 003	Enhances sperm motility	[9]
18	NSC 697923	Reduces sperm motility	This study
19	Pha-665752	Modulates P4 evoked Ca ²⁺ response	[28]
20	Ouabain	Reduces sperm motility	[29]
21	Trequinsin hydrochloride	Increases hyperactivated sperm; increases penetration of sperm into viscous media, increases motility	[28]
22	Leelamine hydrochloride	Modulates P4 evoked Ca ²⁺ response	[28]
23	SLM 6031434 hydrochloride	Induces acrosome reaction	This study
24	JQ1	Interferes with spermatogenesis; reduces sperm count; reduces motility	[30]
25	n-[4-(4-Chlorophenyl)-1,3-thiazol-2-yl] guanidine	Related to Hv1	[27]
26	Nafamostat	Interferes with liquefaction; reduces sperm motility	[31, 32]
27	Anandamide	Reduces sperm motility	[26]
28	Torin 2	Enhances sperm motility	[9]
29	PD 146176	Antagonizing negative effects of ROS on sperm function such as acrosome reaction	[33]
30	Mibefradil	Blocks progesterone induced Ca ²⁺ response	[34]
31	CYM 5442 hydrchloride	Reduces sperm motility; induces acrosome reaction	This study
32	MJ33	Prevents capacitation-associated tyrosine phosphorylation by PKA	[35]
33	HC-056456	Blocks Ca ²⁺ and Na ⁺ entry; blocks onset of hyperactive motility	[36]
34	Urinastatin	Interferes with liquefaction; reduces sperm motility	[31]
35	Pristimerin	Modulates sperm motility; reduces progesterone induced Ca ²⁺ response	[37, 38]
36	CZC 25146	LRRK2 neuronal toxicity	[39]
37	NNC 55-0396 dihydrochloride	Blocks ccl-20 induced Ca ²⁺ response; Modulates CatSper	[34]
38	666-15	Reduces sperm motility; induces acrosome reaction	This study
39	Phenanthroline	Inhibit proteolytic activity of PSA	[32, 40]
40	CCR6 inhibitor	Patent	[41]
41	Oprea1_544341	Reduces sperm motility; role of glycolysis in sperm motility	[42, 43]
42	A23187	Induces acrosome reaction	[44]
43	LRRK2-IN-1	Enhances	[9]
44	Disulfiram	immotile sperm; spermicidal	[8, 45, 46]
45	Forskolin	Stimulates acrosome reaction	[47]

(Continued)

Table 1. (Continued)

STB Number	Name	Published Phenotype	References
46	Alfuzosin	Interferes with ejaculation	[48]
47	Compound 19	Blocks TSSK; similarity to brigatinib	[49]
48	IBMX	Modulates capacitation	[35]
49	Tamsulosin	Interferes with ejaculation	[48]
50	Imiquimod	Reduction of sperm motility; only sperm with X chromosome	[24]
51	Bithionol	Blocks sAC	[15, 17]
52	GF 109203X	Blocks progesterone induced acrosome reaction; blocks progesterone induced Ca ²⁺ response	[18]
53	CTFRinh-172	Blocks progesterone induced acrosome reaction; blocks sperm hyperactivation; blocks recombinant ZP3 induced acrosome reaction	[50]
54	Sodium orthovanadate	Blocks progesterone induced acrosome reaction; blocks progesterone induced Ca ²⁺ response	[18]
55	FPL 64176	Modulates evoked Ca ²⁺ response	[28]
56	183	Reduces acrosome reaction	[51]
57	Esomeprazole	Reduces sperm motility	[52]
58	68853496	Patent; sAC inhibitor	[53]
59	Gossypol	Causes azoospermia	[54]
60	Clofarabine	Enhances sperm motility	[9]
61	B07 hydrochloride	Reduction of sperm motility	[55]
62	KT 5720	Evokes Ca ²⁺ signaling in sperm	[56]
63	Linsitinib	enhances sperm motility	[9]
64	Tafenoquine	Reduces sperm motility; induces acrosome reaction	This study
65	Hesperadin hydrochloride	Reduces sperm motility	This study
66	UK 78282 hydrochloride	Modulates Ca ²⁺ response	[28]
67	Tacrolimus	Blocks progesterone induced acrosome reaction	[18]
68	BPO-27	CTFR inhibitor patent	[57]
69	Tak-063	Enhances sperm motility	[9]
70	LRE1	Prevents hyperactivation	[58]
71	TAE226	Structural similarity with brigatinib (reduces motility)	This study
72	GW 843682X	Enhances sperm motility	[9]
73	GSK1904529A	Inhibition of hyperactivation	[59]
74	AZD 7762 hydrochloride	Reduces sperm motility	This study
75	Cyclosporin A	Blocks progesterone induced acrosome reaction	[18]
76	NVP-AEW541	Inhibition of hyperactivity	[59]
77	PF 431396	Blocks capacitation	[60]
78	JX 401	Modulates Ca ²⁺ response	[28]
79	H 89 dihydrochloride	Evokes Ca ²⁺ signaling in sperm cells	[56]
80	Ro 106-9920	Reduces sperm motility	This study
81	CCR6 inhibitor 1	Potentially blocks acrosome reaction	[34, 61]
82	Herbimycin a	Blocks progesterone induced acrosome reaction; blocks progesterone induced Ca ²⁺ response	[18]
83	Miglustat	Interferes with spermatogenesis	[62]

<https://doi.org/10.1371/journal.pone.0297666.t001>

To compare the effects of STB compounds on specific sperm functions (e.g. motility, acrosome reaction, viability), we screened them in multiple phenotypic functional assays described below. Each sperm function assay has been repeated in multiple donor pools (see [Materials and methods](#)), where each pool consists of cells donated from at least 3 different donors to address donor to donor variability. Each compound has been tested at two concentrations, 10 and 30 μ M, and normalized to vehicle controls (DMSO, see [Materials and methods](#)). We have

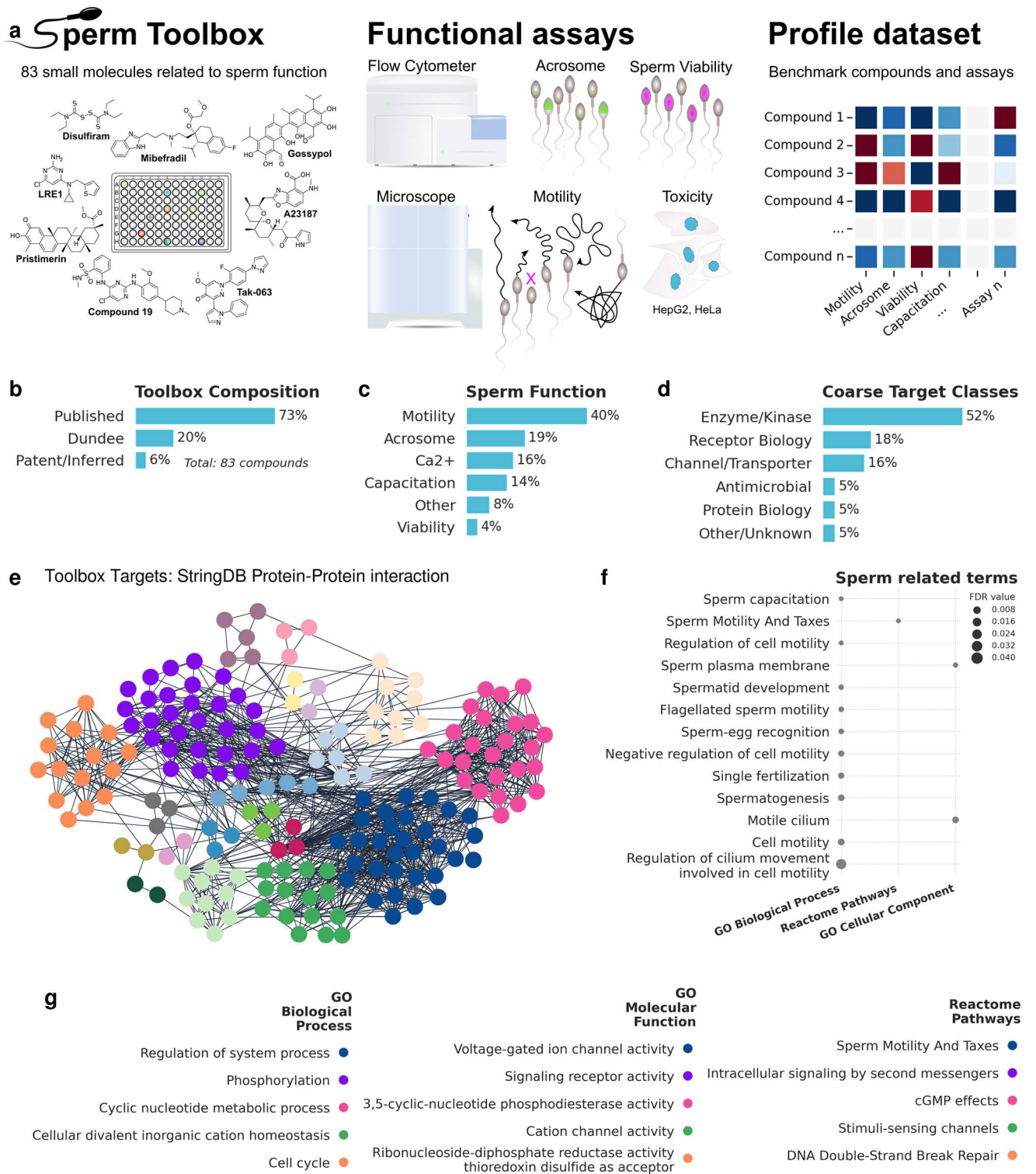


Fig 1. Design and assay of the Sperm Toolbox. (A) The Sperm Toolbox is a collection of 83 small molecules with published, observed or inferred effect on multiple biological processes related to the biogenesis or function of spermatozoa. Human sperm cells have been tested in various assays, which allows comparison of compounds under the same experimental conditions. (B) Summary of the STB compound compositions. (C) Summary of annotated (i.e. published, observed, inferred) sperm function of Toolbox compounds. (D) Coarse classification of STB target classes. (E) StringDB Protein-Protein interaction networks (Full network, confidence 0.7) of putative targets of STB compounds. Protein targets have been inferred from the repurposing hub. Network has been clustered using MCL clustering in Cytoscape. (F) Sperm related terms (GeneOntology, Reactome) within the network shown in (E).

Marker size indicates false discovery rate (FDR). (G) Enriched terms within the biggest subclusters in (E). Redundant terms have been filtered using Cytoscape (cutoff 0.5). Color matches subcluster color in (E).

<https://doi.org/10.1371/journal.pone.0297666.g001>

modified our assays from the previously published workflow [8] by increasing incubation time and introducing conditions that support capacitation, which is a sperm maturation process required for fertilization. These changes expand the range of targets that we can explore in our assays. Our live-sperm-based high-throughput motility assays use an automated microscope to record short timelapse movies. These assays can be run in different protocols (S1 Fig). The first protocol, denoted ‘motility assay’, uses conditions that do not support capacitation and model the state of spermatozoa after ejaculation. In this protocol, cells are incubated with compounds for 10 min prior to the first motility read-out. A second and third motility read-out were performed after 3h and 6h incubation. The second protocol, denoted ‘capacitation assay’, utilises conditions that support capacitation (S1 Fig). In this protocol, cells are kept in conditions supporting capacitation for 30 minutes prior to 30 minutes incubation with compound, followed by motility read-out. Both protocols allow sperm motility to be measured (0.5 sec time-lapse movies allow tracking of cells and calculation of sperm kinematics). In addition, the capacitation assay also scores hyperactivation (i.e. distinct movement patterns concomitant with capacitation which are crucial for fertilization) (Fig 2A).

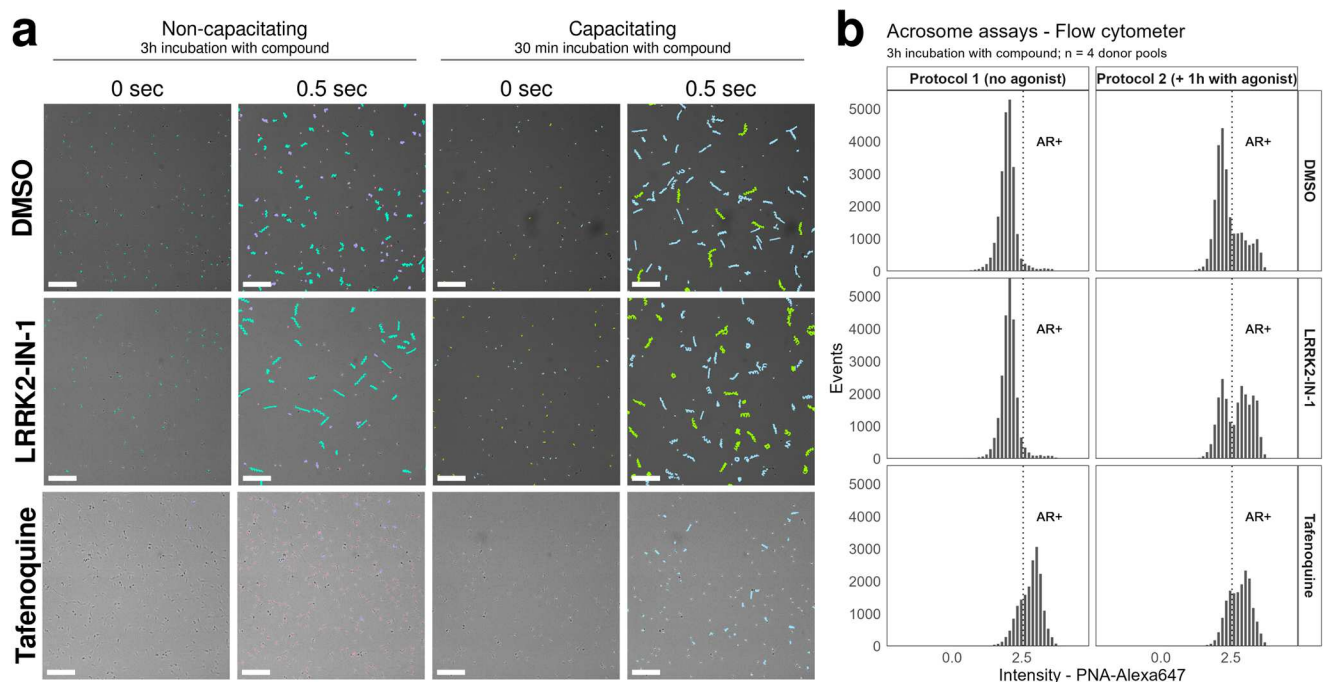


Fig 2. Example assay data from STB compounds. (A) Motility assay example images. Tracked sperm cells treated with DMSO (top row), LRRK2-IN-1 (middle row) or Tafenoquine (bottom row). First (0 sec) and last frame (0.5 sec) of time-lapse movies are shown. Non-capacitating conditions (3h incubation with compound) are shown next to capacitating conditions (30 min incubation with compound). Compound concentration was 10 μ M for all shown conditions. Colors in the non-capacitating assay indicate immotile sperm (pink), non-progressively motile sperm (lavender) or progressively motile sperm (teal). Colors in the capacitation assay indicate hyperactive sperm (green) or non-hyperactive sperm (blue). Scale bars 50 μ m. (B) Acrosome assay examples. Histogram depicting number of events and shift in intensity in the PNA acrosome channel (AR+ population). Acrosome assay protocol 1 (3h incubation with compound, no agonist addition; left panel) and acrosome assay protocol 2 (3h incubation with compound, 1h with agonist; right panel). Comparing DMSO (top row), LRRK2-IN-1 (middle row) and Tafenoquine (bottom row). Compound concentration was 10 μ M for all shown conditions. Dashed line indicates gate (AR+ population).

<https://doi.org/10.1371/journal.pone.0297666.g002>

The acrosome assays utilize a high-throughput flow cytometer (S1 Fig) to score for compounds inducing the acrosome reaction, using peanut agglutinin conjugated to a fluorescent dye to label the outer acrosomal membrane. We have developed two different screening protocols, named ‘acrosome 1’ and ‘acrosome 2’ (S1 Fig). Both protocols are capable of measuring induction of the acrosome reaction upon compound incubation. However, in the acrosome 2 assay, a mixture of physiological agonists of the acrosome reaction is used after cells have been incubated with compounds for 3h under conditions that support capacitation [63]. This addition allowed measurement of inhibition of acrosome reaction (Fig 2B).

We also tested the effect of STB compounds on sperm viability using a complementary live-cell flow-cytometer approach and propidium iodide, which does not permeate intact cell membranes (S1 Fig). In this assay, cells are incubated for 3h with compound in non-capacitating media, prior to viability assessment.

Finally, we also profiled each STB compound for general cytotoxicity in two cell-lines (HeLa and HepG2). The HeLa-based assay uses the cell painting assay [11], while the HepG2 assay utilizes a resazurin based viability read-out [64]. To ensure usability of STB compounds, we also measured aqueous solubility. An overview of additional assays is given in S1 Fig.

This multi-assay compound profiling dataset reveals a broader picture of the action of individual compounds than has been available before and helps find phenotypic similarities induced by compounds. Using unsupervised clustering we observed clusters of compounds with a significant reduction of sperm motility and/or increased amounts of acrosome reacted cells (e.g. Alexidine dihydrochloride, Tafenoquin or A23187) (Fig 3A). Some of these compounds also show significant levels of cytotoxicity in HeLa and/or HepG2 cells as well as decreased sperm viability suggesting that their effect on motility and acrosome reaction is due to toxic effects on spermatozoa. However, other clusters show a significantly increased acrosome reaction and hyperactive motility but do not show effects on sperm viability or general cytotoxicity (e.g. LRRK2-IN-1, Clofarabine, Linsitinib, Trequinsin). One example, LRRK2-IN-1, which enhances motility (progressive and hyperactivated), and increases the number of acrosome reacted cells in our acrosome assay 2 (Fig 2) has been described as a selective inhibitor of LRRK2 kinase, which is involved in Parkinson’s disease [65]. This compound shares some structural similarities with two kinase inhibitors in the STB, which decrease sperm motility: Compound 19, described as a potent Testis-Specific Serine/Threonine Kinase 2 (TSSK2) inhibitor [49] and CZC 25146, annotated as a potent LRRK2 inhibitor (Fig 3B and 3C). These three compounds show different activity profiles in our profiling assays, suggesting they may target different pathways in sperm. Interestingly, LRRK2-IN-1 only increases acrosome reacted cells in conjunction with the agonist cocktail (acrosome assay 2), while without there is no change, highlighting the potential of our assays to pick up different effects of the same compounds, when the assay condition has been slightly modified (Fig 2B).

Interestingly, we observed that a number of compounds show weak or no effects in our assays (S2 Fig; grey datapoints). Given these compounds were sourced from publications, patents, and screening systems that mostly report activity in sperm, we have retained them in the toolbox as they may score in other assays or at other concentrations. Alternatively if they do not score in assays of sperm function, their inclusion serves to correct previous results. Either way, this highlights a key strength of our approach. The STB, and the use of high-throughput assays allows the profiling of compounds using the same assay conditions (i.e. buffer system, incubation time and screening concentration) on a large number of human spermatozoa, thus providing reproducible, quantitative assay results. This combination also facilitates the comparison of the effects of the STB compounds between different physiological states (S1 Fig), in a standardised, rapid, and quantitative manner that has not been previously possible in the field of sperm biology. This standardized approach provides

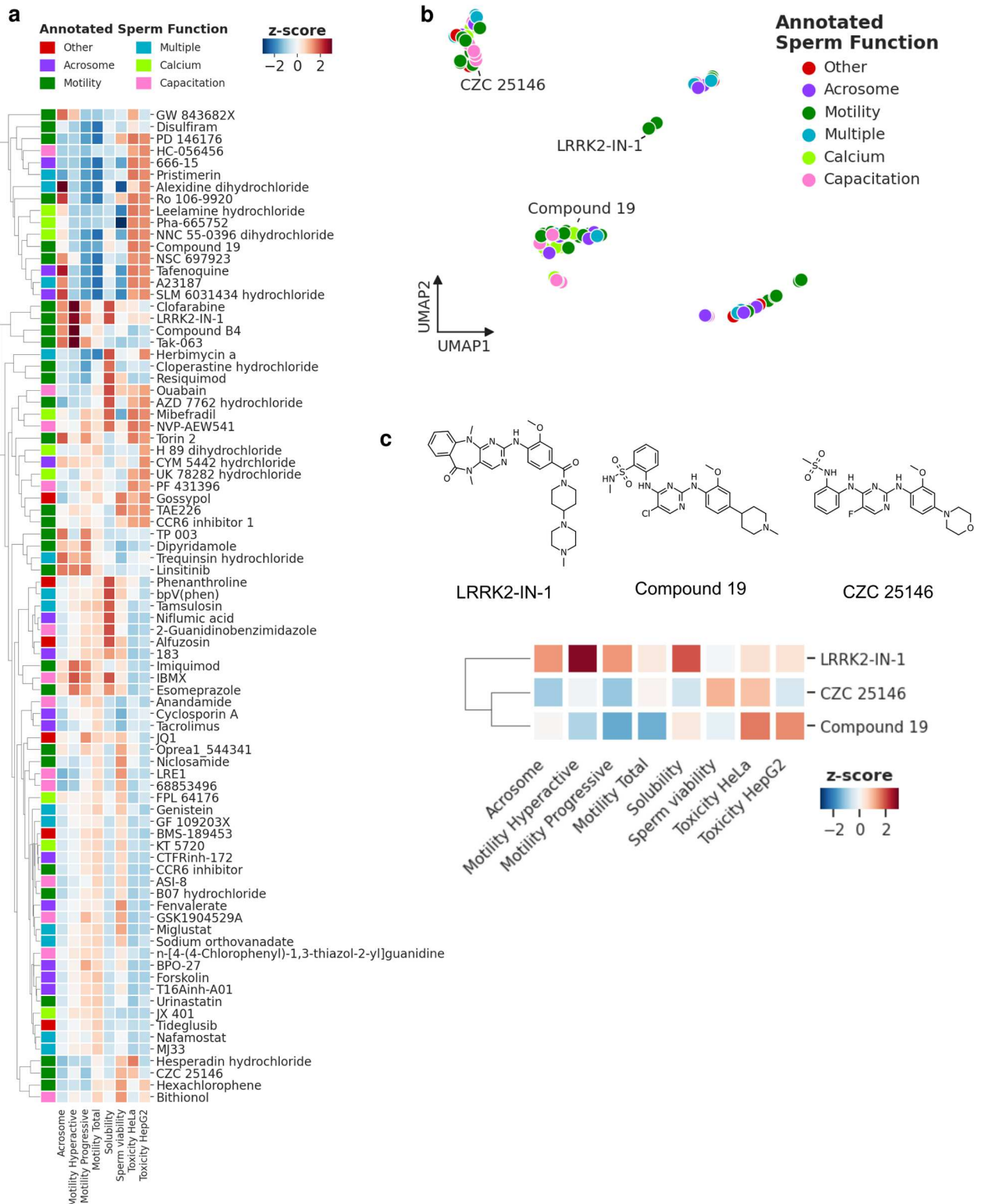


Fig 3. Visualization of assay results on STB compounds. (A) Heatmap showing standardized results (z-score) for every Toolbox compound screened at 30 μ M in selected assays: Acrosome 2 (3h with compound then induction with acrosome reaction using an agonist mix), Motility Hyperactive (30 min capacitation then 30 min incubation with compound), Motility Progressive (3h incubation with compounds under non-capacitating conditions), Motility Total (3h incubation with compound under non-capacitating conditions)), Solubility (aqueous solubility of compounds after 24h incubation), Sperm viability (3h with compound then propidium iodide assay), Toxicity HeLa (24h compound incubation

then cell painting assay), Toxicity HepG2 (70h incubation with compound then resaruzin based assay). Heatmap color indicates decrease (blue) or increase (red). (B) UMAP plot showing groups of compounds with similar assay profile. Color indicates annotated sperm function. (C) Structures and zoom-in of three compounds with related structures described to modulate LRRK2 kinase but appear to have different effects on sperm function.

<https://doi.org/10.1371/journal.pone.0297666.g003>

an excellent starting point to establish control compounds and assay conditions for novel drug discovery campaigns related to contraception and infertility. We provide the complete dataset (see [Supporting information](#)) as well as web interface for data analysis (<https://spermttoolbox.shinyapps.io/spermttoolbox/>).

Conclusion

This work represents the first version of the STB. We aim to release annual updates with newly published compounds that modulate the function of spermatozoa. In our next release we aim to include compounds that were published during the collection of the data for this first release e.g., a sAC inhibitor TDI-11861 [4], a CatSper inhibitor RU1968 [66], a SLO3 inhibitor [67] or other contributions from the research community. The STB is a reference for compounds that modulate sperm function and can become an important resource for benchmarking new assays and new compounds. The STB allows comparison of compounds side-by-side using the same assay conditions therefore mitigating the effects of different conditions e.g., buffers, cell populations which may give erroneous results. We aim to use the STB as a reference in our screens and to define as many factors targeted by STB compounds as possible that are critical for sperm function.

Materials and methods

Ethical approval

Written consent was obtained from each donor in accordance with under local ethical approval (University of Dundee, SMED REC 20/45). Donors for the study were sourced from July 1, 2020 until 31 March, 2023.

Sperm handling

Donated human spermatozoa with normal semen characteristics (concentration, motility; WHO 2010) from healthy volunteers with no known fertility problems were used in this study. Samples were obtained by masturbation after ≥ 48 h of sexual abstinence, liquefied for 30 min at 37°C, and processed using density gradient centrifugation, to separate cells using 40/80% Percoll (Sigma Aldrich, UK) fractions. Every experiment has been performed on pooled spermatozoa from ≥ 3 donors each.

For experiments in non-capacitating conditions, we used a buffer system slightly modified from [23] using Minimal Essential Medium Eagle (Sigma), supplemented with HEPES (1 M solution, Gibco), sodium lactate (DL-solution, Sigma), sodium pyruvate (100 mM solution, Gibco) and bovine serum albumin (7.5% solution, Sigma).

For experiments in capacitating conditions, we used a buffer system slightly modified from the HTF from [63] to a final composition of 97.8 mM NaCl, 4.69 mM KCl, 0.2 mM MgSO₄, 0.37 mM KH₂PO₄, 2.04 mM CaCl₂, 0.33 mM Na-pyruvate, 12.4 mM Na-lactate, 2.78 mM glucose, 21 mM HEPES, 25 mM NaHCO₃ and 3 mg/mL BSA. All buffer components were supplied by Sigma-Aldrich unless otherwise stated and all buffer systems were adjusted to pH 7.4.

Compound handling

Compounds were purchased from vendors indicated in the source publications. 10 mM stock solutions in DMSO were generated. Quality control of each compound has been performed (HPLC-MS and NMR).

Compound plates were generated using 10 mM stock solutions of each STB compound in Echo 384-LDV plates (LP-0200). Prior to running experiments, assay ready plates have been generated using Echo acoustic dispensers (550/555). Each compound was screened at 10 and 30 μ M final concentration. Plate types used for assays: Motility, HeLa cell painting and Sperm viability (PerkinElmer 384-well PhenoPlate), Acrosome (Greiner 384-PP V-bottom plates).

Motility assays

Motility assay (Non-capacitating). Cells were rested for 1 h post density gradient centrifugation. 10 μ L assay buffer was dispensed into assay ready plates, plates were put onto a shaker instrument for 30 seconds using 1,000 rpm, followed by a short pulse centrifugation. Plates were incubated for 15–30 mins at 37°C prior to addition of 10 μ L of pooled cells (at 1 M/mL concentration). Plates were then incubated in a Yokogawa CV7000 microscope set at 37°C for 10 min prior to imaging. Imaging a 384-well plate requires \leq 20 min using the following settings: 2 positions per well, 0.5 sec timelapse movies (3x binning; 11 ms exposure time, 22 ms interval time, 500 intervals). A second and a third timepoint was recorded after 3 h, and 6 h incubation at 37°C. Data has been processed as described in [8]. Assay results are from 5 separate donor pools.

Capacitation assay. Cells were capacitated for 30 min at 37°C and 5% CO₂. 10 μ L assay buffer was dispensed into assay ready plates, plates were put onto a shaker instrument for 30 seconds using 1,000 rpm, followed by a short pulse centrifugation. Plates were incubated for 15–30 mins at 37°C and 5% CO₂ prior to addition of 10 μ L of pooled cells (at 1 M/mL concentration). This was followed by a 20 min compound incubation step at 37°C and 5% CO₂. Plates were then incubated in our CV7000 microscope set at 37°C and 5% CO₂ for 10 min prior to imaging. Imaging a 384-well plate requires \leq 20 min using the following settings: 2 positions per well, 0.5 sec timelapse movies (3x binning; 11 ms exposure time, 22 ms interval time, 500 intervals). Data has been processed as described in [8]. Assay results are from 8 separate donor pools.

Acrosome assays. In both experimental protocols, 20 μ L of pooled cells (at 1 M/mL concentration) were capacitated for 3 h in presence of compounds at 37°C and 5% CO₂. In protocol a, cells were then stained with PNA-Alexa647 (ThermoFisher, 1 mg/mL, final dilution 1:1000) and Hoechst dye for 30 min at 37°C and 5% CO₂, followed by 10 min fixation at RT by adding 20 μ L of 4% paraformaldehyde (Sigma) in PBS. In mode b, a cocktail of progesterone, prostaglandin A and NH₄Cl was added to all compound wells and incubated for an additional 1 h at 37°C and 5% CO₂ to induce acrosome reaction. This was followed by fixation and staining protocol as for mode a. In both protocols fixatives and stains were removed and cells were resuspended in PBS, using a 405 platewasher. Plates were then measured on a Novocyte Advanteon high-throughput flow cytometer. Gating of data has been performed using the flow cytometer software to export counts of gated populations. The Acrosome 1 protocol using 3h incubation with compounds was run in 7 separate donor pools. The Acrosome 2 protocol using 3h incubation with compound followed by 1h stimulation of acrosome reaction was run in 5 separate donor pools.

Sperm viability assay. Cells were incubated with compounds for a total of 3h at 37°C. 20 min prior, propidium iodide was added (1.5 mM solution, Sigma, final dilution 1:2000). Plates

were analysed using a Satorius iQue Screener with settings to allow for a sampling time of ≤ 30 min per 384-well plate. Gating of data has been performed using the flow cytometer software to export counts of gated populations. This assay was tested in 5 separate donor pools.

HeLa cell painting. Cells were incubated with compounds for 24h. Cell painting was performed after [11] with modifications suggested in [68]:

- 10 μL of MitoTracker (Life Technologies, M22426, 1 mM stock in DMSO as recommended by the vendor) is added to 50 μL of media (1:2000 final dilution)
- Phalloidin-Alexa594 (Life Technologies, A12381, 66 μM stock in methanol as recommended by the vendor, 1:500 final dilution)
- WGA-Alexa594 (Life Technologies, W11262, 5 mg/mL stock, 1:250 final dilution) was used
- DAPI (5 mg/mL stock; 1:1000 final dilution) was used to stain DNA
- 20 μL of staining solution was used
- PBS was used for washing steps

HepG2 cytotoxicity assay. HepG2 assay was performed as described previously [64].

Aqueous solubility. The aqueous solubility of test compounds was measured using UHPLC, as previously described in [69].

Data analysis. For cell painting imaging data, plates of images were imported into OMERO [70]. Cell segmentation and feature extraction were performed using CellProfiler 3 [71] and using instructions from [11] to normalize data. Timelapse data of motility assays was processed as previously described [8]. Acrosome flow cytometry data was processed using NovoCyte Express (PerkinElmer), to gate out background debris, select for single cells and calculate percentage values of acrosome positive events. Sperm viability flow cytometry data was processed using iQue Forecyt (Satorius) to gate out background debris, select for single cells and calculate percentage values of propidium iodide positive events. For all sperm functional assays, data was normalized to DMSO control wells. Each plate had 16 DMSO wells. A median value of all 16 wells was used to calculate percentage change ($\text{pct_ch} = ((\text{value} / \text{DMSO_median}) - 1) * 100$). For sperm viability we calculated fold-change ($\text{fc} = 1 - (\text{value} / \text{DMSO_median})$). For each compound we calculated a p-value using a Welch's t-test implemented in scipy (`scipy.stats.ttest_ind`) comparing a compound's effect to DMSO controls.

The following packages have been used for data analysis and visualization:

R: Tidyverse, plotly, pHeatmap

Python: Pandas, NumPy, SciPy, Seaborn, Matplotlib, UMAP

Supporting information

S1 Fig. STB assay protocols. (a) Diagram showing sperm motility assays run on high-content microscope with compound incubation time/read-out time. Motility assays were run under non-capacitating conditions, with a 1 hour recovery phase after preparation of spermatozoa. Capacitation assays were run under capacitating conditions. Spermatozoa are capacitated for 30 min prior to compound incubation. (b) Diagram showing assays run on high-throughput flow cytometer to measure acrosome status and sperm viability. Incubation time with compound, agonist and read-out times are indicated. Acrosome assays run under capacitating

conditions. Sperm viability assay run under non-capacitating condition. (c) Additional assays performed on sperm toolbox compounds.

(TIF)

S2 Fig. Volcano plots of acrosome, capacitation and motility assays. Upper panel indicates 10 μ M, lower plot panels indicate 30 μ M. Incubation times for assays are 3h with compound for Acrosome 2 assay and Motility (non-capacitating) (PM = progressive motility), and 30 min for the Capacitation assay (HA = hyperactive motility, PM = progressive motility). Dotted line indicates significance levels of 0.05. Any point above dotted line is indicated in magenta. Buffer conditions are indicated below Plot title (capacitating vs. non-capacitating).

(TIF)

S1 Table. STB assay data. Supporting data underlying the findings in this article.

(CSV)

Acknowledgments

The authors thank Alistair Langlands, John Raynor and the UK National Phenotypic Screening Centre for support and assistance performing the automated screening assays.

Author Contributions

Conceptualization: Kevin D. Read, Christopher L. R. Barratt, Ian H. Gilbert, Jason R. Swedlow.

Data curation: Anthony Richardson, Zoe C. Johnston, Rachel Myles, Neil R. Norcross.

Formal analysis: Franz S. Gruber.

Funding acquisition: Kevin D. Read, Sarah Martins da Silva, Christopher L. R. Barratt, Ian H. Gilbert, Jason R. Swedlow.

Investigation: Franz S. Gruber, Zoe C. Johnston, Neil R. Norcross, David P. Day, Irene Georgiou, Laura Sesma-Sanz, Caroline Wilson, Kevin D. Read.

Methodology: Franz S. Gruber, Zoe C. Johnston, Rachel Myles, Neil R. Norcross, David P. Day, Irene Georgiou, Laura Sesma-Sanz, Caroline Wilson.

Project administration: Christopher L. R. Barratt, Ian H. Gilbert, Jason R. Swedlow.

Software: Franz S. Gruber.

Supervision: Kevin D. Read, Sarah Martins da Silva, Christopher L. R. Barratt, Ian H. Gilbert, Jason R. Swedlow.

Validation: Franz S. Gruber, David P. Day, Irene Georgiou, Laura Sesma-Sanz, Caroline Wilson, Kevin D. Read.

Visualization: Franz S. Gruber, Jason R. Swedlow.

Writing – original draft: Franz S. Gruber, Zoe C. Johnston, Jason R. Swedlow.

Writing – review & editing: Franz S. Gruber, Anthony Richardson, Zoe C. Johnston, Kevin D. Read, Sarah Martins da Silva, Christopher L. R. Barratt, Ian H. Gilbert, Jason R. Swedlow.

References

1. Norcross NR, Georgiou I, Johnston ZC, Gruber FS, Swedlow JR, Read KD, et al. Male contraceptive development: A medicinal chemistry perspective. *Eur J Med Chem.* 2022; 243:114709. <https://doi.org/10.1016/j.ejmech.2022.114709> PMID: 36087385
2. Anderson DJ, Johnston DS. A brief history and future prospects of contraception. *Science.* 2023; 380(6641):154–8. <https://doi.org/10.1126/science.adf9341> PMID: 37053322
3. Service CA, Puri D, Hsieh TC, Patel DP. Emerging concepts in male contraception: a narrative review of novel, hormonal and non-hormonal options. *Ther Adv Reprod Health.* 2023; 17:26334941221138323. <https://doi.org/10.1177/26334941221138323> PMID: 36909934
4. Balbach M, Rossetti T, Ferreira J, Ghanem L, Ritagliati C, Myers RW, et al. On-demand male contraception via acute inhibition of soluble adenylyl cyclase. *Nat Commun.* 2023; 14(1):637. <https://doi.org/10.1038/s41467-023-36119-6> PMID: 36788210
5. Faber EB, Wang N, John K, Sun L, Wong HL, Burban D, et al. Screening through Lead Optimization of High Affinity, Allosteric Cyclin-Dependent Kinase 2 (CDK2) Inhibitors as Male Contraceptives That Reduce Sperm Counts in Mice. *J Med Chem.* 2023; 66(3):1928–40. <https://doi.org/10.1021/acs.jmedchem.2c01731> PMID: 36701569
6. Carlson EJ, Francis R, Liu Y, Li P, Lyon M, Santi CM, et al. Discovery and Characterization of Multiple Classes of Human CatSper Blockers. *ChemMedChem.* 2022; 17(15):e202000499. <https://doi.org/10.1002/cmdc.202000499> PMID: 35644882
7. Hau RK, Tash JS, Georg GI, Wright SH, Cherrington NJ. Physiological Characterization of the Transporter-Mediated Uptake of the Reversible Male Contraceptive H2-Gamendazole Across the Blood-Testis Barrier. *J Pharmacol Exp Ther.* 2022; 382(3):299–312. <https://doi.org/10.1124/jpet.122.001195> PMID: 35779861
8. Gruber FS, Johnston ZC, Barratt CL, Andrews PD. A phenotypic screening platform utilising human spermatozoa identifies compounds with contraceptive activity. *Elife.* 2020; 9. <https://doi.org/10.7554/eLife.51739> PMID: 31987071
9. Gruber FS, Johnston ZC, Norcross NR, Georgiou I, Wilson C, Read KD, et al. Compounds enhancing human sperm motility identified using a high-throughput phenotypic screening platform. *Hum Reprod.* 2022; 37(3):466–75. <https://doi.org/10.1093/humrep/deac007> PMID: 35048946
10. Perlman ZE, Slack MD, Feng Y, Mitchison TJ, Wu LF, Altschuler SJ. Multidimensional Drug Profiling By Automated Microscopy. *Science.* 2004; 306(5699):1194–8. <https://doi.org/10.1126/science.1100709> PMID: 15539606
11. Bray MA, Singh S, Han H, Davis CT, Borgeson B, Hartland C, et al. Cell Painting, a high-content image-based assay for morphological profiling using multiplexed fluorescent dyes. *Nat Protoc.* 2016; 11(9):1757–74. <https://doi.org/10.1038/nprot.2016.105> PMID: 27560178
12. Samby K, Besson D, Dutta A, Patra B, Doy A, Glossop P, et al. The Pandemic Response Box horizontal line Accelerating Drug Discovery Efforts after Disease Outbreaks. *ACS Infect Dis.* 2022; 8(4):713–20.
13. Brown PJ, Muller S. Open access chemical probes for epigenetic targets. *Future Med Chem.* 2015; 7(14):1901–17. <https://doi.org/10.4155/fmc.15.127> PMID: 26397018
14. Corsello SM, Bittker JA, Liu Z, Gould J, McCarren P, Hirschman JE, et al. The Drug Repurposing Hub: a next-generation drug library and information resource. *Nat Med.* 2017; 23(4):405–8. <https://doi.org/10.1038/nm.4306> PMID: 28388612
15. Kleinboelting S, Ramos-Espiritu L, Buck H, Colis L, van den Heuvel J, Glickman JF, et al. Bithionol Potently Inhibits Human Soluble Adenylyl Cyclase through Binding to the Allosteric Activator Site. *J Biol Chem.* 2016; 291(18):9776–84. <https://doi.org/10.1074/jbc.M115.708255> PMID: 26961873
16. Saalau-Bethell SM, Berdini V, Cleasby A, Congreve M, Coyle JE, Lock V, et al. Crystal structure of human soluble adenylyl cyclase reveals a distinct, highly flexible allosteric bicarbonate binding pocket. *ChemMedChem.* 2014; 9(4):823–32. <https://doi.org/10.1002/cmdc.201300480> PMID: 24616449
17. Balbach M, Fushimi M, Huggins DJ, Steegborn C, Meinke PT, Levin LR, et al. Optimization of lead compounds into on-demand, nonhormonal contraceptives: leveraging a public-private drug discovery institute collaboration dagger. *Biol Reprod.* 2020; 103(2):176–82.
18. Baron L, Fara K, Zapata-Carmona H, Zuniga L, Kong M, Signorelli J, et al. Participation of protein kinases and phosphatases in the progesterone-induced acrosome reaction and calcium influx in human spermatozoa. *Andrology.* 2016; 4(6):1073–83.
19. Chen Z, Shu N, Wang Y, Yang Y, Shao Z, Tian F, et al. Tideglusib, a prospective alternative to nonoxynol-9 contraceptive. *Contracept X.* 2019; 1:100007. <https://doi.org/10.1016/j.conx.2019.100007> PMID: 32550526

20. Orta G, Ferreira G, Jose O, Trevino CL, Beltran C, Darszon A. Human spermatozoa possess a calcium-dependent chloride channel that may participate in the acrosomal reaction. *J Physiol*. 2012; 590(11):2659–75. <https://doi.org/10.1113/jphysiol.2011.224485> PMID: 22473777
21. Vega SG, Guzman P, Garcia L, Espinosa J, Cortinas de Nava C. Sperm shape abnormality and urine mutagenicity in mice treated with niclosamide. *Mutat Res*. 1988; 204(2):269–76. [https://doi.org/10.1016/0165-1218\(88\)90099-7](https://doi.org/10.1016/0165-1218(88)90099-7) PMID: 3278217
22. Chung SS, Wang X, Roberts SS, Griffey SM, Reczek PR, Wolgemuth DJ. Oral administration of a retinoic Acid receptor antagonist reversibly inhibits spermatogenesis in mice. *Endocrinology*. 2011; 152(6):2492–502. <https://doi.org/10.1210/en.2010-0941> PMID: 21505053
23. Tardif S, Madamidola OA, Brown SG, Frame L, Lefievre L, Wyatt PG, et al. Clinically relevant enhancement of human sperm motility using compounds with reported phosphodiesterase inhibitor activity. *Hum Reprod*. 2014; 29(10):2123–35. <https://doi.org/10.1093/humrep/deu196> PMID: 25124668
24. Umehara T, Tsujita N, Shimada M. Activation of Toll-like receptor 7/8 encoded by the X chromosome alters sperm motility and provides a novel simple technology for sexing sperm. *PLoS Biol*. 2019; 17(8): e3000398. <https://doi.org/10.1371/journal.pbio.3000398> PMID: 31408454
25. O'Rand MG, Silva EJ, Hamil KG. Non-hormonal male contraception: A review and development of an Eppin based contraceptive. *Pharmacol Ther*. 2016; 157:105–11. <https://doi.org/10.1016/j.pharmthera.2015.11.004> PMID: 26593445
26. Mishra AK, Kumar A, Yadav S, Anand M, Yadav B, Nigam R, et al. Functional insights into voltage gated proton channel (Hv1) in bull spermatozoa. *Theriogenology*. 2019; 136:118–30. <https://doi.org/10.1016/j.theriogenology.2019.06.015> PMID: 31254725
27. Hong L, Singh V, Wulff H, Tombola F. Interrogation of the intersubunit interface of the open Hv1 proton channel with a probe of allosteric coupling. *Sci Rep*. 2015; 5:14077. <https://doi.org/10.1038/srep14077> PMID: 26365828
28. McBrinn RC, Fraser J, Hope AG, Gray DW, Barratt CLR, Martins da Silva SJ, et al. Novel pharmacological actions of trequinsin hydrochloride improve human sperm cell motility and function. *Br J Pharmacol*. 2019; 176(23):4521–36. <https://doi.org/10.1111/bph.14814> PMID: 31368510
29. Syeda SS, Sanchez G, Hong KH, Hawkinson JE, Georg GI, Blanco G. Design, Synthesis, and in Vitro and in Vivo Evaluation of Ouabain Analogues as Potent and Selective Na,K-ATPase alpha4 Isoform Inhibitors for Male Contraception. *J Med Chem*. 2018; 61(5):1800–20.
30. Matzuk MM, McKeown MR, Filippakopoulos P, Li Q, Ma L, Agno JE, et al. Small-molecule inhibition of BRDT for male contraception. *Cell*. 2012; 150(4):673–84. <https://doi.org/10.1016/j.cell.2012.06.045> PMID: 22901802
31. Matsuda Y, Oshio S, Yazaki T, Umeda T, Akihama S. The effect of some proteinase inhibitors on liquefaction of human semen. *Hum Reprod*. 1994; 9(4):664–8. <https://doi.org/10.1093/oxfordjournals.humrep.a138568> PMID: 8046020
32. Anamthathmakula P, Winuthayanon W. Mechanism of semen liquefaction and its potential for a novel non-hormonal contraception dagger. *Biol Reprod*. 2020; 103(2):411–26.
33. Walters JLH, De Iulius GN, Dun MD, Aitken RJ, McLaughlin EA, Nixon B, et al. Pharmacological inhibition of arachidonate 15-lipoxygenase protects human spermatozoa against oxidative stress. *Biol Reprod*. 2018; 98(6):784–94. <https://doi.org/10.1093/biolre/iy058> PMID: 29546268
34. Diao R, Wang T, Fok KL, Li X, Ruan Y, Yu MK, et al. CCR6 is required for ligand-induced CatSper activation in human sperm. *Oncotarget*. 2017; 8(53):91445–58. <https://doi.org/10.18632/oncotarget.20651> PMID: 29207656
35. Lee D, Moawad AR, Morielli T, Fernandez MC, O'Flaherty C. Peroxiredoxins prevent oxidative stress during human sperm capacitation. *Mol Hum Reprod*. 2017; 23(2):106–15. <https://doi.org/10.1093/molehr/gaw081> PMID: 28025393
36. Carlson AE, Burnett LA, del Camino D, Quill TA, Hille B, Chong JA, et al. Pharmacological targeting of native CatSper channels reveals a required role in maintenance of sperm hyperactivation. *PLoS One*. 2009; 4(8):e6844. <https://doi.org/10.1371/journal.pone.0006844> PMID: 19718436
37. Mannowetz N, Miller MR, Lishko PV. Regulation of the sperm calcium channel CatSper by endogenous steroids and plant triterpenoids. *Proc Natl Acad Sci U S A*. 2017; 114(22):5743–8. <https://doi.org/10.1073/pnas.1700367114> PMID: 28507119
38. Brenker C, Schiffer C, Wagner IV, Tuttmann F, Ropke A, Rennhack A, et al. Action of steroids and plant triterpenoids on CatSper Ca(2+) channels in human sperm. *Proc Natl Acad Sci U S A*. 2018; 115(3):E344–E6. <https://doi.org/10.1073/pnas.1717929115> PMID: 29305558
39. Ramsden N, Perrin J, Ren Z, Lee BD, Zinn N, Dawson VL, et al. Chemoproteomics-based design of potent LRRK2-selective lead compounds that attenuate Parkinson's disease-related toxicity in human neurons. *ACS Chem Biol*. 2011; 6(10):1021–8. <https://doi.org/10.1021/cb2002413> PMID: 21812418

40. Robert M, Gibbs BF, Jacobson E, Gagnon C. Characterization of prostate-specific antigen proteolytic activity on its major physiological substrate, the sperm motility inhibitor precursor/seminalogelin I. *Biochemistry*. 1997; 36(13):3811–9. <https://doi.org/10.1021/bi9626158> PMID: 9092810
41. Dairaghi D DD, Kalisiak J, Lange CW, Leleti MR, Li Y, Lui RM, et al. ChemoCentryx Inc, assignee. CCR6 compounds. United States. US 20190167659A1 2018. 11 p. <https://patents.google.com/patent/US20190167659A1/en>.
42. Danshina PV, Qu W, Temple BR, Rojas RJ, Miley MJ, Machius M, et al. Structural analyses to identify selective inhibitors of glyceraldehyde 3-phosphate dehydrogenase-S, a sperm-specific glycolytic enzyme. *Mol Hum Reprod*. 2016; 22(6):410–26. <https://doi.org/10.1093/molehr/gaw016> PMID: 26921398
43. Miki K, Qu W, Goulding EH, Willis WD, Bunch DO, Strader LF, et al. Glyceraldehyde 3-phosphate dehydrogenase-S, a sperm-specific glycolytic enzyme, is required for sperm motility and male fertility. *Proc Natl Acad Sci U S A*. 2004; 101(47):16501–6. <https://doi.org/10.1073/pnas.0407708101> PMID: 15546993
44. Aitken RJ, Buckingham DW, Fang HG. Analysis of the responses of human spermatozoa to A23187 employing a novel technique for assessing the acrosome reaction. *J Androl*. 1993; 14(2):132–41. PMID: 8514619
45. Lal N, Jangir S, Bala V, Mandalapu D, Sarswat A, Kumar L, et al. Role of disulfide linkage in action of bis (dialkylaminethiocarbonyl)disulfides as potent double-Edged microbicidal spermicide: Design, synthesis and biology. *Eur J Med Chem*. 2016; 115:275–90. <https://doi.org/10.1016/j.ejmech.2016.03.012> PMID: 27084496
46. Danscher G, Rebbe H. Effects of two chelating agents, oxine and diethyldithiocarbamate (Antabuse), on stainability and motility of human sperms. *J Histochem Cytochem*. 1974; 22(10):981–5. <https://doi.org/10.1177/22.10.981> PMID: 4139184
47. Doherty CM, Tarchala SM, Radwanska E, De Jonge CJ. Characterization of two second messenger pathways and their interactions in eliciting the human sperm acrosome reaction. *J Androl*. 1995; 16(1):36–46. PMID: 7768751
48. Hellstrom WJ, Sikka SC. Effects of alfuzosin and tamsulosin on sperm parameters in healthy men: results of a short-term, randomized, double-blind, placebo-controlled, crossover study. *J Androl*. 2009; 30(4):469–74. <https://doi.org/10.2164/jandrol.108.006874> PMID: 19201696
49. Hawkinson JE, Sinville R, Mudaliar D, Shetty J, Ward T, Herr JC, et al. Potent Pyrimidine and Pyrrolo-pyrimidine Inhibitors of Testis-Specific Serine/Threonine Kinase 2 (TSSK2). *ChemMedChem*. 2017; 12(22):1857–65.
50. Li CY, Jiang LY, Chen WY, Li K, Sheng HQ, Ni Y, et al. CFTR is essential for sperm fertilizing capacity and is correlated with sperm quality in humans. *Hum Reprod*. 2010; 25(2):317–27. <https://doi.org/10.1093/humrep/dep406> PMID: 19923167
51. Baggelaar MP, den Dulk H, Florea BI, Fazio D, Perruzza D, Bernabo N, et al. ABHD2 Inhibitor Identified by Activity-Based Protein Profiling Reduces Acrosome Reaction. *ACS Chem Biol*. 2019; 14(12):2943. <https://doi.org/10.1021/acscchembio.9b00824> PMID: 31697462
52. Kumar A, Kumar R, Flanagan J, Langstrom B, Bjorndahl L, Darreh-Shori T. Esomeprazole reduces sperm motility index by targeting the spermic cholinergic machinery: A mechanistic study for the association between use of proton pump inhibitors and reduced sperm motility index. *Biochem Pharmacol*. 2020; 182:114212. <https://doi.org/10.1016/j.bcp.2020.114212> PMID: 32866455
53. Buchmann B HM, Cancho-Grande Y, Kosemund D, Schirok H, Nguyen D, Fritsch M, et al. Bayer Pharma AG, assignee. Azaindoles as inhibitors of soluble adenylylase cyclase. United States. US20090098189A1.2008 Sep 1. 45 p. <https://patents.google.com/patent/US20090098189A1/en>.
54. Coutinho EM. Gossypol: a contraceptive for men. *Contraception*. 2002; 65(4):259–63. [https://doi.org/10.1016/s0010-7824\(02\)00294-9](https://doi.org/10.1016/s0010-7824(02)00294-9) PMID: 12020773
55. Yang M, Zhi R, Lu L, Dong M, Wang Y, Tian F, et al. A CCR5 antagonist-based HIV entry inhibitor exhibited potent spermicidal activity: Potential application for contraception and prevention of HIV sexual transmission. *Eur J Pharm Sci*. 2018; 117:313–20. <https://doi.org/10.1016/j.ejps.2018.02.026> PMID: 29496533
56. Wang T, Young S, Krenz H, Tuttmann F, Ropke A, Krallmann C, et al. The Ca(2+) channel CatSper is not activated by cAMP/PKA signaling but directly affected by chemicals used to probe the action of cAMP and PKA. *J Biol Chem*. 2020; 295(38):13181–93. <https://doi.org/10.1074/jbc.RA120.013218> PMID: 32703901
57. Verkman AS SD, inventor; University of California, assignee. Pyrimido-pyrrolo-oxazine-dione compound inhibitors of the cystic fibrosis transmembrane conductance regulator protein and uses therefor. United States. US20140080821A1.2012 May 25. 10 p. <https://patents.google.com/patent/US20140080821A1/en>.

58. Ramos-Espiritu L, Kleinboelting S, Navarrete FA, Alvau A, Visconti PE, Valsecchi F, et al. Discovery of LRE1 as a specific and allosteric inhibitor of soluble adenylyl cyclase. *Nat Chem Biol.* 2016; 12(10):838–44. <https://doi.org/10.1038/nchembio.2151> PMID: 27547922
59. Wang J, Qi L, Huang S, Zhou T, Guo Y, Wang G, et al. Quantitative phosphoproteomics analysis reveals a key role of insulin growth factor 1 receptor (IGF1R) tyrosine kinase in human sperm capacitation. *Mol Cell Proteomics.* 2015; 14(4):1104–12. <https://doi.org/10.1074/mcp.M114.045468> PMID: 25693802
60. Alvau A, Battistone MA, Gervasi MG, Navarrete FA, Xu X, Sanchez-Cardenas C, et al. The tyrosine kinase FER is responsible for the capacitation-associated increase in tyrosine phosphorylation in murine sperm. *Development.* 2016; 143(13):2325–33. <https://doi.org/10.1242/dev.136499> PMID: 27226326
61. Duan YG, Wehry UP, Buhren BA, Schrupf H, Olah P, Bunemann E, et al. CCL20-CCR6 axis directs sperm-oocyte interaction and its dysregulation correlates/associates with male infertility. *Biol Reprod.* 2020; 103(3):630–42.
62. van der Spoel AC, Jeyakumar M, Butters TD, Charlton HM, Moore HD, Dwek RA, et al. Reversible infertility in male mice after oral administration of alkylated imino sugars: a nonhormonal approach to male contraception. *Proc Natl Acad Sci U S A.* 2002; 99(26):17173–8. <https://doi.org/10.1073/pnas.262586099> PMID: 12477936
63. Prajapati P, Kane S, McBrinn RC, Dean MS, Martins da Silva SJ, Brown SG. Elevated and Sustained Intracellular Calcium Signalling Is Necessary for Efficacious Induction of the Human Sperm Acrosome Reaction. *Int J Mol Sci.* 2022; 23(19). <https://doi.org/10.3390/ijms231911253> PMID: 36232560
64. Baragana B, Hallyburton I, Lee MC, Norcross NR, Grimaldi R, Otto TD, et al. A novel multiple-stage antimalarial agent that inhibits protein synthesis. *Nature.* 2015; 522(7556):315–20. <https://doi.org/10.1038/nature14451> PMID: 26085270
65. Deng X, Dzamko N, Prescott A, Davies P, Liu Q, Yang Q, et al. Characterization of a selective inhibitor of the Parkinson's disease kinase LRRK2. *Nat Chem Biol.* 2011; 7(4):203–5. <https://doi.org/10.1038/nchembio.538> PMID: 21378983
66. Schierling T, Tosi B, Eisenhardt C, Reining S, Daniliuc CG, Brenker C, et al. Synthesis and Functional Characterization of Novel RU1968-Derived CatSper Inhibitors with Reduced Stereochemical Complexity. *ACS Pharmacol Transl Sci.* 2023; 6(1):115–27. <https://doi.org/10.1021/acsptsci.2c00188> PMID: 36654752
67. Lyon M, Li P, Ferreira JJ, Lazarenko RM, Kharade SV, Kramer M, et al. A selective inhibitor of the sperm-specific potassium channel SLO3 impairs human sperm function. *Proc Natl Acad Sci U S A.* 2023; 120(4):e2212338120. <https://doi.org/10.1073/pnas.2212338120> PMID: 36649421
68. Cimini BA, Chandrasekaran SN, Kost-Alimova M, Miller L, Goodale A, Fritchman B, et al. Optimizing the Cell Painting assay for image-based profiling. *Nat Protoc.* 2023; 18(7):1981–2013. <https://doi.org/10.1038/s41596-023-00840-9> PMID: 37344608
69. Smith A, Wall RJ, Patterson S, Rowan T, Rico Vidal E, Stojanovski L, et al. Repositioning of a Diaminopyrimidopyrimidinone Series Confirmed to Target the Cyclin-Dependent Kinase CRK12 for Use in the Treatment of African Animal Trypanosomiasis. *J Med Chem.* 2022; 65(7):5606–24. <https://doi.org/10.1021/acs.jmedchem.1c02104> PMID: 35303411
70. Allan C, Burel JM, Moore J, Blackburn C, Linkert M, Loynton S, et al. OMEERO: flexible, model-driven data management for experimental biology. *Nat Methods.* 2012; 9(3):245–53. <https://doi.org/10.1038/nmeth.1896> PMID: 22373911
71. McQuin C, Goodman A, Chernyshev V, Kamensky L, Cimini BA, Karhohs KW, et al. CellProfiler 3.0: Next-generation image processing for biology. *PLoS Biol.* 2018; 16(7):e2005970. <https://doi.org/10.1371/journal.pbio.2005970> PMID: 29969450